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ENGLISH

You have purchased a RHINO product, which guarantees quality and durability, since the BABA-180M balance was manufactured according to American quality standards. Keep this instruction manual, it will help you understand the operation of the RHINO BABA-180M balance will be useful for future reference.

TECHNICAL SPECIFICATIONS

CAPACITY

400 lb

DIVISION

0.2 lb

POWER

2 pcs x AAA
batteries (included)

TEMPERATURE

32 °F – 104 °F

QUICK GUIDE TO USE

- Remove the scale from its packaging, make sure it is in perfect condition and includes all its accessories (manual and 2 batteries).
- Place the scale on a firm, level surface. Any movement, no matter how small, can affect the accuracy.
- Preferably place your feet in the center of the scale.

PRECAUTIONS

- The display will show **"Lo"** when the battery needs to be replaced.
- The display will show **"Err"** when the scale is overloaded. Remove excess weight to prevent damage.
- Put the scale in a cool and dry place.
- For family use only. DO NOT USE for commercial purposes.
- Clean the product with a soft, slightly damp (not wet) cloth. Do not immerse the scale in water.

WARNINGS

Do not open the scale by yourself, or the warranty will be void.

To clean the scale, use a damp cloth or sponge, never a direct stream of water.

Protect the scale from insects, do not use spray, only solid products against insects.

Do not drop the scale suddenly.

Do not direct airflow from fans or heaters directly at the scale.



Do not expose your scale to dripping or splashing liquids.

INSTALLING THE BATTERY

- Open the cover at the bottom of the scale as shown in the diagram.
- Install the batteries taking into account the polarity indicated in the compartment.
- Replace the battery compartment cover.

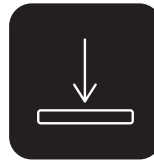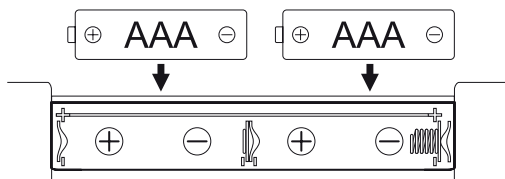
FUNCTIONING

To ensure the accuracy of the scale, place it on a hard, flat surface.

Note. Avoid placing on a slippery surface.


- Carefully step on the scale, it will turn on automatically. Stay on the scale without moving until the weight appears on the screen, when the numbers on the display flash 3 times and immediately remain fixed, it will indicate your weight.

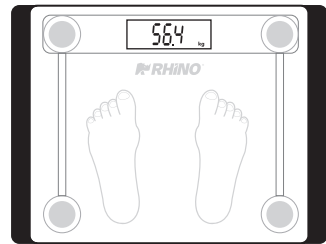UNIT CONVERSION

The unit change button is located on the back of the scale (UNIT), as shown in the diagram. With the (UNIT) button you can change the unit of measure, from kilograms (kg) to pounds (lb) and vice versa.

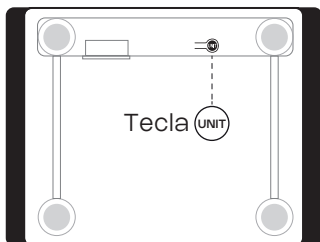WARNING MESSAGE

-  Indicates that there is overweight on the scale.
-  Indicates that the battery is low.
-  If the display shows the letter **"C"** after being weighed, it means that the weight shown on the screen has some error and must be weighed again.

USE WARNING

-  The glass of the scale can be very slippery when wet.
-  Stand tall while weighing yourself on the scale.
-  Do not hit or drop the scale.
-  Handle the scale with care as it may damage the accuracy of the scale.
-  Clean the scale with a damp cloth but avoid getting water inside the scale.

-  Do not use chemicals to clean it.
-  Keep the scale on a cool, ventilated surface.
-  Always keep the scale in a horizontal position.
-  If the scale does not turn on, check that the batteries are installed or the batteries may be low. Replace batteries.
-  If an error appears on the scale or the scale does not turn on for a while, remove the batteries for 5 seconds and then put them back. If the problem continues, contact your dealer.
-  Do not use this scale for commercial purposes.

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WARRANTY

Rhino Maquinaria warrants to customers who purchase any Rhino products from an authorized seller that Rhino products will be free of defects in material or workmanship for the duration of 2 years from the date of purchase.

During the 2-years warranty period, Rhino Maquinaria will replace the product that has a manufacturing defect as long as it has been determined by Rhino Maquinaria.

This warranty does not apply to defects caused by misuse or alterations of the product, nor does it apply if the damage was caused by repairs or service provided or attempted by anyone other than Rhino Maquinaria.

To make this warranty effective, you must first contact our Customer Service by dialing (619) 393-1079 or by mail to info@rhinoglobal.net, once you have contacted our Customer Service they will request from you a warranty claim by email.

Rhino Maquinaria will review the warranty claim. If the warranty claim is found to be valid, please return the product with a copy of this invoice and order ID to the office of invoice and order ID to the indicate address.

Upon receiving the product, Rhino Maquinaria will send a replacement product to customer's address.

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Rhino Maquinaria S.A. de C.V.